



Synthesis Report

Tracer : FDG Duo c (1.0.0.1)

Customer batch number

Date : 07 Apr 2023

Operator : Administrator (Administrator)
Cassette test sequence : GE duoFDGc PETCS003FL

GUID: 5009A3BE-3871-4BD4-8437-70807D88A1E9

Synthesis sequence : GE duoFDGc LTV PETCS020FL B

GUID: BF8E77DE-BF5C-41C3-BEF2-F3F657BCC330

Synthesizer ID : Module FL2-0288 Software version : 3.3.0.14 - PLC: 7.75

Workflow ID: : 794-2

Calibration report number : Id n°1034 FDG Duo c (1.0.0.1) 22 Jun 2021

12:49:29

Cassette expiry date : 01 Jul 2023

Cassette reference : 030004062307231052

Process report

Synthesizer test : Successful.
Cassette test : Successful.
Synthesis : Successful.
Pause : Not Used.
Manual Mode : Not Used.

Cassette rinsing : Please refer to rinsing report 794.

The values of radioactivity detectors are only indicative.

Yield summary

Time at reagent connection : 07 Apr 2023 04:00:27
Time at start of synthesis : 07 Apr 2023 09:13:11
Time at end of synthesis : 07 Apr 2023 09:35:26

Synthesis duration : 22 min. 15 s.

Target volume : 5.0 ml

Starting activity : 401567 MBq I Tracer activity : 285002 MBq

I at end of synthesis

Tracer activity : N/A

at end of transfer

Volume of tracer : 22 ml

Concentration : 12955 MBg/ml

: 71 %

: 81 %

end of synthesis


I Estimated yield

I at end of synthesis

Estimated yield

decay corrected



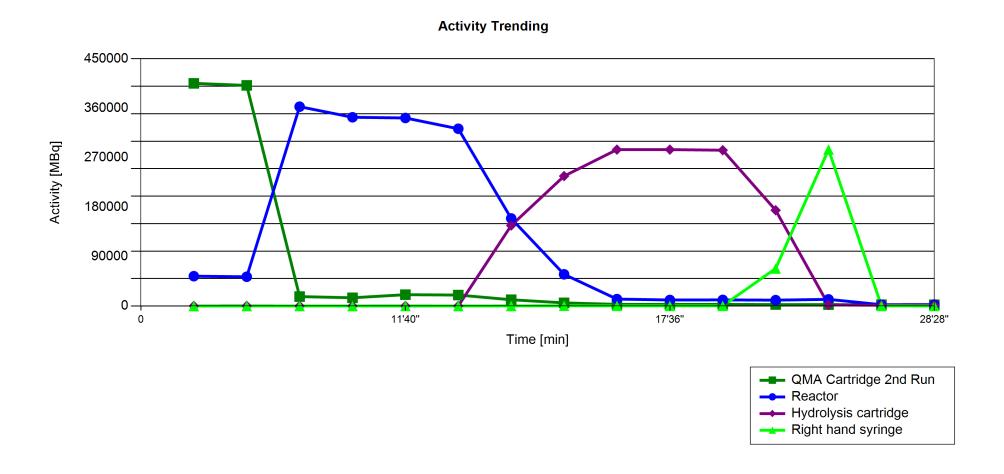



Trending

Measurments Performed	Time from start	QMA Cartridge 2nd Run [MBq]	Reactor [MBq]	Hydrolysis cartridge [MBq]	Right hand syringe [MBq]	Heater reactor [°C]
FDG run2 - Start of Synthesis	00'02"	405353	54308	120	0	25
FDG run2 - 18F trapped on QMA	00'49"	401567	53265	281	53	25
FDG run2 - 18F eluted from QMA to reactor	06'53"	17182	362886	120	53	120
FDG run2 - 18F dried	11'19"	15041	343579	80	161	125
FDG run2 - Start of labelling	11'40"	20655	342317	161	269	124
FDG run2 - End of labelling	13'40"	20002	322695	161	107	125
FDG run2 - First trapping on tC18	15'00"	11437	159431	146579	215	105
FDG run2 - Second trapping on tC18	15'51"	5744	57703	236540	511	97
FDG run2 - Third trapping on tC18	16'47"	3029	12612	284742	403	91
FDG run2 - tC18 cartridge washed	17'36"	3055	11011	284742	349	87
FDG run2 - Start of basic hydrolysis	18'05"	3029	11157	283615	376	85
FDG run2 - End of basic hydrolysis	20'25"	2689	10745	174405	68304	76
FDG run2 - FDG transferred to Syringe S3	22'15"	2559	12152	2013	285002	70
FDG run2 - Transferred to Product Collection Vial	28'22"	2376	2449	1449	295	56
Cassette rinsing	28'28"	2219	2643	1610	107	56

The values of radioactivity detectors are only indicative.

Log

Severity	Event Time	Message
Information	07 Apr 2023 00:48:21	Cassette RFID tag successfully read. (#1173)
Information	07 Apr 2023 05:19:54	Synthesis succeeded #1



Information

07 Apr 2023 09:42:55

Synthesis succeeded #2



Operator signature :

Responsible signature :